

4050X Gynaecological Examination Table

Hygienic easy-to-use examination table for multiple purposes



Main features:

1. Easy to adjust entirely hands-free, which makes healthcare personnel's work more hygienic and effective
2. Designed especially for gynaecological examinations, but it's suitable also for minor procedures and general examinations
3. New features that improve usability, comfort and patient safety

The 4050X gynaecological examination table is the most versatile examination table in our range. It is particularly well suited for gynaecological examinations, but it can also be used in minor procedures and general examinations. The 4050X gynaecological examination table features a three-piece table top with seat section, middle section and back section adjustments. Seat, middle section and height can be adjusted from the front of the table using the new gynaecological foot control set, which makes healthcare personnel's work more hygienic and effective. As an option, you can also add a two-piece hands-free adjustment bar to the side of the table or use the hand control that is a standard feature.

In addition to the new hygienic gynaecological foot control set, the upholstery material is extremely durable and antimicrobial fabric that is also easy to clean, waterproof and has excellent disinfectant features. As standard the 4050X is equipped with adjustable heel supports (stepless and no levers) and a gynaecological bowl, that is easy to move out of the way when not needed. 4050X's motors lower using gravity and in emergency situations - if there's something caught between the mechanism – table rises automatically.

There are two models of 4050X available. The 4050X M has a smaller number of features, whereas 4050X F is more extensive - f.ex. it has a hand control with memory slot, electric seat adjustment and the new gynaecological foot control set (see technical data and features for more details).

LOJER® For easy care


Technical data and features	4050X F	4050X M
Total width	75 cm	75 cm
Patient surface width	70 cm	70 cm
Total length	200 cm	200 cm
Total weight	120 kg	120 kg
Height adjustment	49-97 cm	49-97 cm
Back section length	70 cm	70 cm
Back section adjustment	0-73°, electrical	0-73°, electrical
Middle section length	82 cm	82 cm
Middle section adjustment	0-73°, electrical	0-73°, gas spring
Seat section length	39 cm	39 cm
Seat section adjustment	0-25°	Non-adjustable
Electrical Trendelenburg / anti-Trendelenburg (optional)	20° / -10°	20° / -10°
Castors	Ø 125 mm locking castors	Ø 125 mm locking castors
Other standard features	Hand control with memory slot and simultaneous movement of different sections gynaecological foot control set for height, seat and middle section adjustments steplessly adjustable heel supports gynaecological bowl accessory rails in seat section paper roll stand and cutter safety feature: angle between seat and middle section min. 105° safety feature: auto-lifting in emergency situations safety feature: lock function by pressing two pedals on the gynaecological foot control set or buttons on hand control for 3 sec	Hand control with safety lock function by pressing two buttons for 3 sec steplessly adjustable heel supports gynaecological bowl accessory rails in seat section paper roll stand and cutter safety feature: auto-lifting in emergency situations
Upholstery	Several material and colour options, more information on request	Several material and colour options, more information on request
Safe working load (SWL)	210 kg	210 kg
Electric system	24V	24V
CE -marking	Yes	Yes
Warranty	2 years	2 years
Country of Manufacture	Finland	Finland

Options

- Trendelenburg and anti-Trendelenburg incl. battery
- Central locking castors Ø125 or Ø150 mm
- Central locking conductive twin castors Ø 150 mm
- Two-piece hands-free adjustment bar for height and back section
- Foot pedal for height and back section adjustment
- Battery powered

Accessories

- Pair of calf supports, white plastic (for minor procedures). Includes rapid clamps and fastening straps.
- Pair of calf supports, black moulded cushion (for long-term procedures). Includes rotation clamps and fastening straps.
- IV-pole straight, bended or antimicrobial coating
- Accessory rail
- Accessory rail, brass
- Rapid clamp Ø 20 mm
- Side rail
- Side rail, antimicrobial coating
- Anaesthesia frame
- Push handles
- Examination lamp (attachment to accessory rail or to the frame)
- Neck cushion
- Arm support for sampling
- Adapter for fitting an IV-stand (4 possible positions)

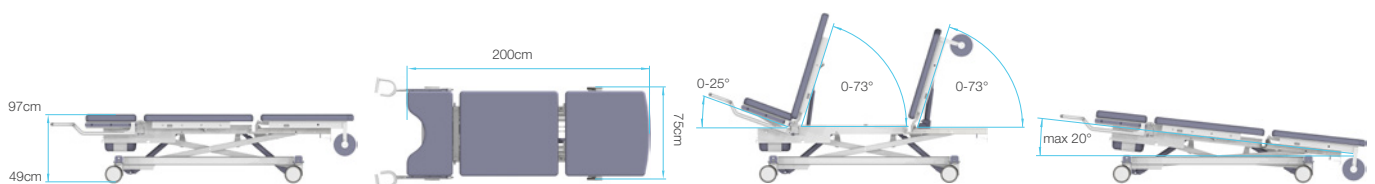

LOJER® For easy care